

An integrated surgical solution

Stille and Dr. Charles Randquist team-up for a state-of-the-art breast augmentation solution.

When one of highest quality surgical instruments teams-up with one of the most innovative breast augmentation surgeons in the world, you can only expect a winning outcome. That's exactly what Stille, in cooperation with Dr. Charles Randquist in Stockholm, Sweden have done.

When properly implemented, the procedure and instrument set can allow for shorter operating times and increased patient turnover. Couple that with fewer

complications, less bleeding and faster recovery and you have a better solution for your breast augmentation practise.

Efficiency In Numbers

The Stille-Randquist partnership has not only created a more well thought, stream lined surgical procedure but also a more efficient set of surgical instruments.



| Stille-Randquist
Instruments |

Stille - Surgical Perfection

Stille is one of the oldest medical device companies in the world. We develop, manufacture and sell proprietary high-quality surgical instruments and specialty tables to leading surgeons worldwide. We have done so for more than 170 years.

Stille's offering consists of a strong brand and products renowned for their superior quality. We manage our business in two distinct business areas—Surgery and Patient Positioning.

Stille surgical instruments are considered state of the art by many surgeons around the world. All Stille instruments are handcrafted as an art by highly skilled instrument makers in Eskilstuna, Sweden.

The slender design and excellent balance are the characteristics which make the Stille Instruments easily recognizable. Their highly polished surface and delicate design make Stille instruments lighter and easier to keep clean, thus more durable against discoloration and corrosion. All this provides Stille instruments with a unique and outstanding feeling.

Art. no: 950-99101-1 2EN03 © 2012-10-08 Stille. All rights reserved.



Stille-Randquist breast augmentation set

Efficiency and increased patient turnover



Surgical Perfection. For life.

Stille AB

Gustav III:s Boulevard 42
SE-169 73 Solna, Sweden
Phone: +46 8 588 580 00
Fax: +46 8 588 580 05
info@stille.se

www.stille.se

Stille Surgical Inc.

943 Parkview Blvd
Lombard, IL 60148
United States
Phone: 1.800.6651614

Distributed by:

VERKAUF **TapMed** 
Medizintechnik Handels GmbH
Gewerbepark 10
34317 Halbachswald-Ehlen
GERMANY
Tel. +49 (0)5606 53060-0
Fax +49 (0)5606 53060-44
www.tapmed.de



Surgical Perfection. For life.

Bringing breast augmentation to perfection



The Stille-Randquist breast augmentation set

Description	Article code
STILLE-RANDQUIST Complete instruments set	199-99101
RULER gold plated 15 cm 6 in	199-30010-154
SCALPEL HANDLE no. 3 std, black coated handle	199-50010-303
VOLKMANN-STILLE RETRACTOR 4 sharp prongs, black coated, 22 cm 8.75 in	108-24680-403
STILLE BREAST RETRACTOR double-ended, black coated handle, 23 cm 9 in	108-26038
BJÖRK-STILLE DIATHERMY FORCEPS insulated, 25 cm 10 in	102-9845-25
OLSEN-HEGAR NEEDLEHOLDER with scissors, black coated, TC-inserts 18 cm 7 in	184-15590-183
STILLE-RAPP KNOT-HOLDING FORCEPS, black coated, 13 cm 5 in	104-15140-3
MAYO-STILLE SUPERCUT DISSECTING SCISSORS, curved, 15 cm 6 in	101-8083-15
STILLE TISSUE FORCEPS, black coated, 1x2 teeth, 15 cm 6 in	102-9060-153
CRILE-STILLE ARTERY FORCEPS, black coated, curved, 16 cm 6.25 in	103-10182-163
INSTRUMENT FOLDER, three fold	224-04-0215

The Randquist procedure and the Stille instrument set can allow for shorter operating times and increased patient turnover

When used in tandem with the Randquist surgical procedure, the innovative Stille-Randquist set of surgical tools represent the next level for breast augmentation by delivering more predictable results, thus achieving a more positive patient experience.

An integrated Surgical Solution


By letting a better surgical procedure influence the design of the instrument set, we've been able to provide the industry with a more completely integrated surgical solution that goes beyond what is currently considered state-of-the-art.

- Put simply, the Randquist procedure is an efficient solution that includes:
- ✓ Enhanced, simplified preoperative planning eliminating the need for sizers
 - ✓ Creating a non-traumatic, sharp coagulating dissection thus eliminating all of the perforators and minimizing bleeding.
 - ✓ Performing all surgical tasks from one side of the table
 - ✓ One-touch glove change

Randquist Procedure Overview


The efficiency of the Randquist procedure not only supports higher-level breast augmentation procedures, but it can also contribute to shortening the operative time.

When implemented properly, as outlined in the "Procedure Overview" below, it can lead to procedures that can be performed within 20 minutes, contributing to a more positive patient experience.

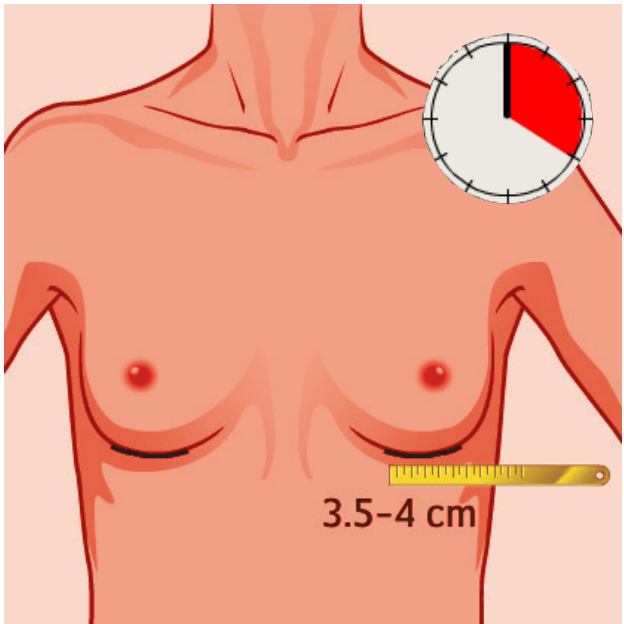

Procedure Overview step 1-6

- 1 Re-measure and mark the pre-marked surgical area with the **Golden ruler**.
- 2 Make the incision through first tissue layer with the **Scalpel**. Hold the wound open during initial dissection with the **Sharp Volkmann retractor**.
- 3 Dissect second tissue layer down to fascia with a dissection needle. Hold the wound open during further pocket creation with the **Double sided breast retractor**.

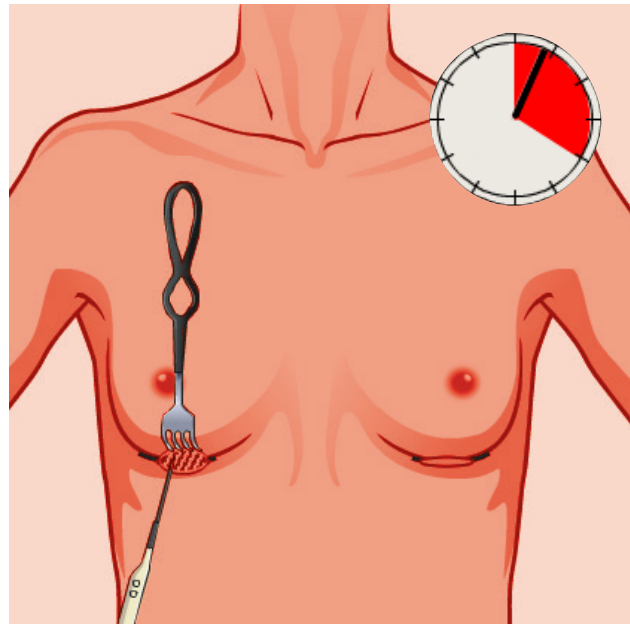

- 4 Create a pocket for the implant with monopolar diathermy. Use diathermy on possible bleedings during dissection with the **Insulated Stille dressing forceps**
- 5 Suture with the **Needle-holder scissors**.
- 6 Hold the suture knot when suturing the lower tissue layer with the **Stille Rapp Knot holding forceps**.

